



PANTHERA™ FUSION BAR



Engineered shape, form and function based on dental science.



700+
platforms

Unmatched implant
compatibility

25
YEAR

Warranty
on the implant bar

Bar Benefits

30% STRONGER

Cantilever tests have shown that the fracture resistance of the Zirconia supported by a Titanium Implant Bar was 30% higher compared to unsupported Zirconia¹.

LOW PROFILE BAR

Little to no bone reduction is required.

ROUNDED CONTOURS AT THE TOP OF THE BAR

No sharp angles to avoid pressure points.

HYGIENIC DESIGN

Precisely designed to replicate your patient's planned intaglio.

OVERLAY MATERIAL - THE CHOICE IS YOURS!

Zirconia or PMMA, multi-layer or monolithic. Select the material that best suits your patients' needs.



WWW.PANTHERADENTAL.COM

Please do not hesitate to contact our support team
if you need any assistance.

INFO@PANTHERADENTAL.COM | 1.418.527.0388




1. These tests were performed internally by Panthera Dental using industry standard protocols and methodology.




PANTHERA™ FUSION BAR



Two innovative bar profiles designed to support
all types of fixed prosthesis.

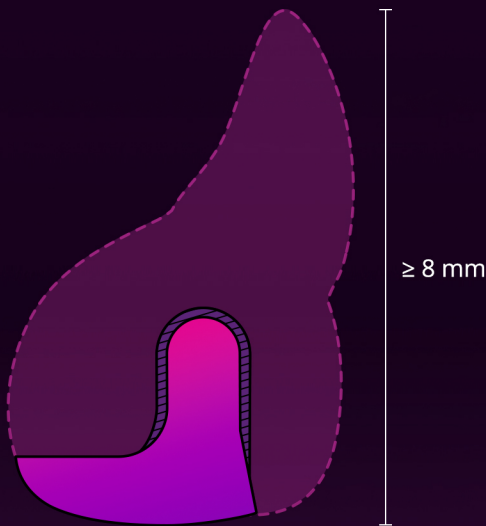

PFB™

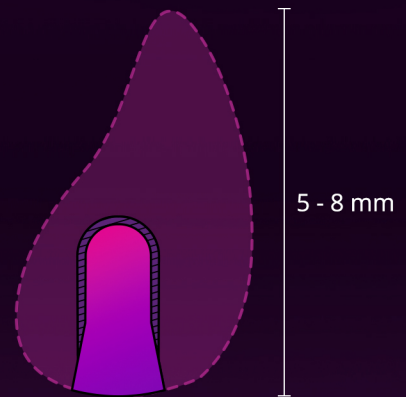

PFB™ Bridge

| FP2-FP3 min. 5 mm | Prosthesis type Low profile bar | FP1-FP2 min. 3.6 mm |
|----------------------|--|---------------------------|
| Buccal | 10° chamfer | Dual : buccal and lingual |
| ✓ | 90° lingual seating | |
| ✓ | Constant 100 micron cement gap | ✓ |
| ✓ | Vertical seating exact positioning | ✓ |
| ✓ | Intimate fit of the overlay on the bar | ✓ |
| ✓ | Low transition line | ✓ |
| ✓ | Angulated screw channels | ✓ |



WWW.PANTHERADENTAL.COM

Please do not hesitate to contact our support team
if you need any assistance.

INFO@PANTHERADENTAL.COM | 1.418.527.0388






PANTHERA™ FUSION BAR



WORK YOUR WAY! Quick turnaround - Easy as 1 2 3.

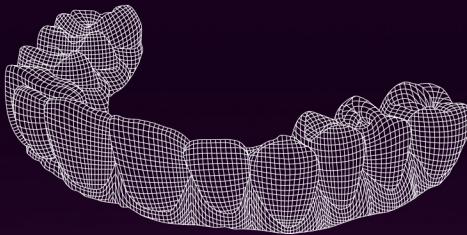

SYNCHRONIZED DIGITAL WORKFLOW

- 1** Panthera splits and creates the bar design in the try-in / prototype and sends you the STL file of the overlay.
- 2** While Panthera mills the bar, you can simultaneously mill the overlay in zirconia or any other material.
- 3** You receive the bar, you cement it with your characterized overlay and voilà!



BAR + ZIRCONIA MILLING SERVICE

- 1** Panthera splits and creates the bar design in the try-in / prototype.
- 2** Panthera mills both the bar and the overlay in zirconia.
- 3** You receive the bar and the overlay, you characterize it, cement it and voilà!



WWW.PANTHERADENTAL.COM

Please do not hesitate to contact our support team
if you need any assistance.

INFO@PANTHERADENTAL.COM | 1.418.527.0388





PANTHERA™ FUSION BAR



Synchronized Digital Workflow

Introducing the first splitting protocol that allows in-lab milling of the overlay while Panthera mills the implant bar.



WWW.PANTHERADENTAL.COM

Please do not hesitate to contact our support team
if you need any assistance.

INFO@PANTHERADENTAL.COM | 1.418.527.0388

